Brushy

# SMART KIDS TOOTHBRUSH

AGU SKT6

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EN USER MANUAL

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## CONTENTS

1	Introduction.....	31
2	List of symbols .....	32
3	Scope of application.....	32
4	Complete set.....	32
5	Basic specifications.....	33
6	Device description .....	34
7	Safety instructions.....	34
8	Before you start.....	36
9	Operating principles and procedures .....	37
10	Device cleaning .....	39
11	Troubleshooting.....	39
12	Storage, transportation and operation rules.....	40
13	Disposal.....	40
14	Certification.....	40
15	Manufacturer's warranty.....	40

## 1 INTRODUCTION

Dear Customer! Thank you for having chosen our products!

Smart kids toothbrush **AGU SKT6** is an electric toothbrush designed specifically for effective cleaning of children's teeth.

To use all the toothbrush functions, install the AGU Brushy app on the phone that your child will be using. You can also get statistics and reports on your child's teeth brushing on another device by installing the AGU.

With Smart kids toothbrush **AGU SKT6**, children can easily learn how to brush their teeth properly. Smart kids toothbrush has a bright design and an interactive training application, making the tooth brushing an exciting game for your child. Using the AGU Brushy app, your child will learn how to brush teeth properly and will develop a healthy habit in a short time. And parents will be able to see the statistics and recommendations for brushing teeth for their child.

The replacement head is designed with the children's teeth structure features in mind. Soft bristle provides extremely safe brushing of children's teeth and gums. Smart kids toothbrush **AGU SKT6** is recommended for use by children aged 3 and older. The brushing procedure must be always carried out under the supervision of adults.

### Smart kids toothbrush AGU SKT6 key functions:

- training interactive application for Android/iOS;
- built-in gyroscope and accelerometer for accurate tracking of the brush head location and movement in the mouth;
- built-in gyroscope and accelerometer for tracking brush movement;
- charging the device via USB-C port;
- IPX6 water ingress protection.

### Oral area brushing timer (Smart Brushy Timer)

The standard brushing time is 2 minutes. Every 20 seconds the toothbrush stops vibrating for 0.5 seconds, indicating that the brushing area should be changed.

### Easy start function

The intensity of vibration gradually increases during the first few brushing cycles to provide comfort when using a new toothbrush.

Please read this instruction manual carefully before using the device.

## 2 LIST OF SYMBOLS

Symbol      Meaning

- The product complies with the basic requirements of:
- Radio Equipment Directive 2014/53/EU;
  - Directive 2011/65/EU RoHS on the restriction of the use of certain hazardous substances in electrical and electronic equipment.

WEEE (Waste Electrical and Electronic Equipment Directive). The symbol on the product or its package means that this product does not fall under the category of domestic waste. To avoid possible harm to the environment and human health, separate such wastes from others and dispose of them in accordance with accepted standards.

- The product complies with the basic requirements of:
- TR CU 020/2011 «On electromagnetic compatibility of technical devices».
  - TR CU 007/2011 «On the safety of products for children and adolescents».

Warning/Attention.

Read the instruction manual before use.

Manufacturer.

## 3 SCOPE OF APPLICATION

Smart kids toothbrush **AGU SKT6** is designed for oral care. The toothbrush is a personal hygiene product and must be used individually. The device is intended for home use. Do not use the device for purposes other than intended.

### WARNING

The toothbrush is not a toy and is not intended for play. Smart kids toothbrush **AGU SKT6** contains small parts which can be swallowed by children. The brushing procedure must be always carried out under the supervision of adults.

## 4 COMPLETE SET

Main unit – 1 pc.

USB-C cable – 1 pc.


Bristle replacement brush head – 2 pcs.

Charging port silicone plug – 1 pc.

Stickers – 1 pc.

User manual – 1 pc.

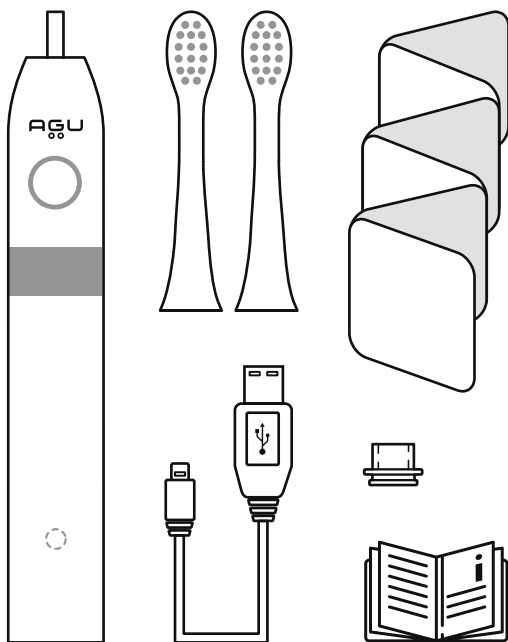
Replacement brush heads for Smart kids toothbrush **AGU SKT6** can be purchased separately.

## 5 BASIC SPECIFICATIONS

Type	Electric toothbrush
Model	AGU SKT6
Power source	Battery ICR10440 320 mAh 3.7 V
Dimensions (L x W x H), mm	180 x 20 x 20
Weight, g	45
BT version	4.0
Maximum vibration frequency (movements per minute)	~11000
Charging current	Up to 250 mA (DC 5 V)
Storage conditions	From +10°C to +30°C with relative humidity up to 65%
Transport conditions	From -10°C to +40°C with relative humidity up to 75%
Operation conditions	From +10°C to +40°C with relative humidity up to 65%
Stiffness bristles	Soft
Service life	2 years

Materials	Replaceable brush head: ABS plastic; main unit: polycarbonate; silicone insert: thermoplastic elastomer; button: polycarbonate; plug: silicone
-----------	---

## 6 DEVICE DESCRIPTION
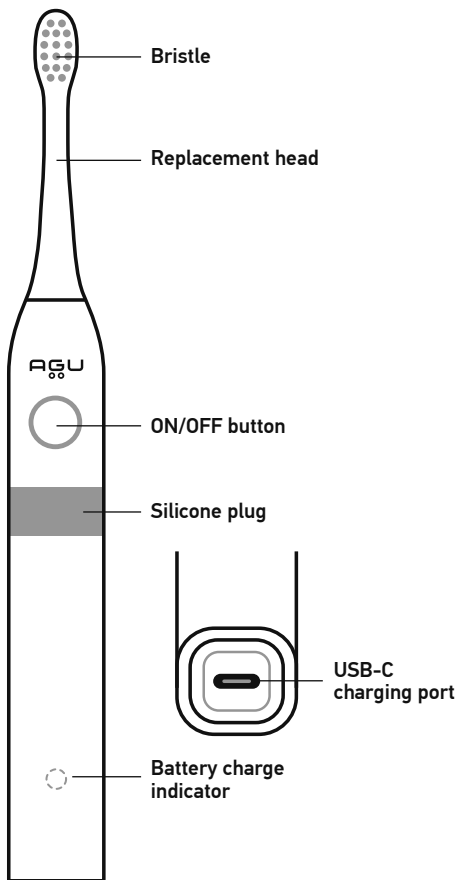
## 7 SAFETY INSTRUCTIONS

### ATTENTION

- The toothbrush is not a toy and is not intended for play. Smart kids toothbrush **AGU SKT6** contains small parts which can be swallowed by children. The brushing procedure must be always carried out under the supervision of adults.
- Basic precautions must always be observed when using any electrical device.

Read the user manual carefully before using the device, and keep it for future reference.

- This device is intended only for household use. Do not use the device for purposes other than intended. Use the device only for the purposes described in the instruction manual. The use of the device for purposes other than those specified in this manual is considered as a violation of the operation rules.
- Periodically check the USB-C cable for possible damage.
- Switch off the device after each use.
- Consult your dentist before using a toothbrush if any of the following concerns your child:
  - sensitive teeth or gums;
  - recent periodontal or dental surgery;
  - undergoing treatment for oral problems.
- This smart kids toothbrush is designed exclusively for oral care in the home. Do not allow children to chew or bite the toothbrush. Do not allow children to play with the toothbrush.
- The toothbrush is a personal hygiene product and must be used individually.

Keep the USB-C cable and power source away from water.

- Do not use a toothbrush while the USB-C cable is connected.
- This device is not suitable for use by people (including children) with impaired physical and mental abilities, and the same people with a lack of experience and knowledge about using this device until they will be instructed about using the device.

### **WARNING**

- Do not lay the USB-C cable near sources of heat, gas or electric cookers. Make sure that the cable does not twist or bend, does not come in contact with sharp objects. Damage to the cable may cause malfunctions not covered by the warranty.
- Never use a damaged cable to charge the device.
- Do not disassemble the device.
- If any malfunctions are detected, stop using the device and contact a specialized service center. Do not attempt to repair the device yourself. All repairs must be carried out by an authorized service center. Unprofessional repairs may result in damage to the device, injury, and property damage.
- The manufacturer shall be not liable for any damage or loss caused by non-observance of safety instructions, improper or unintended use of the device.

Do not use the smart kids toothbrush while taking a bath or shower.

- Stop using the toothbrush and consult your dentist if there is an increased bleeding of the gums or bleeding of the gums after 1 week of use. Visit your dentist if you experience any discomfort or pain while using the toothbrush.

- The smart kids toothbrush head should be replaced every 3 months or more frequently when the bristles become worn out.
- Stop using the device if any part (brush head, handle) is damaged. There are no user-serviceable parts in the device. If the device is damaged, contact a service center.


## 8 BEFORE YOU START

 Read the instruction manual carefully before using the device, and keep it for future reference.


- 1 Get the device and its accessories out of the box, remove all packing materials, promotional stickers and labels.

Your toothbrush handle is waterproof (IPX6), electrically safe and designed specifically for use in the bathroom.


- 2 Connect the USB-C cable to the 5 V / 1 A power adapter (power adapter is not included).

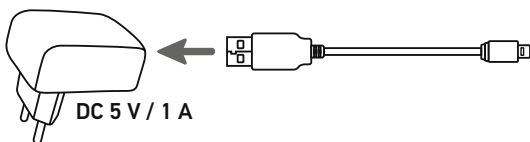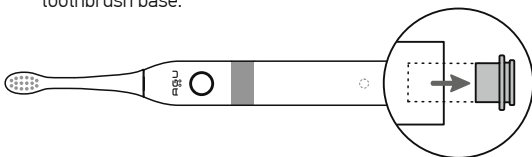- 3 Remove the silicone plug of the charging port from the toothbrush base.

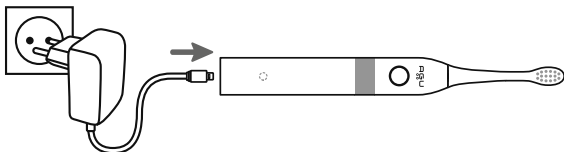- 4 Connect the main unit of the toothbrush to the charging unit.

The charging indicator blinks , while the main unit of the toothbrush is charging. When the main unit is fully charged, the charging indicator goes out . It takes up to 2 hours to fully charge the battery.

- 5 When charging is complete, disconnect the cable and replace the silicone plug. The device is ready to use.

Make sure that the device is completely dry before recharging.

Use only replacement heads recommended by the manufacturer.

## 9 OPERATING PRINCIPLES AND PROCEDURES


### OPERATION WITHOUT USING MOBILE APPLICATION

Wet the brush and squeeze some toothpaste on it.

We start brushing teeth with the upper jaw on the right side.

During teeth brushing we go through three zones consistently (right/middle/left), each in accordance with it's own brushing technique (see pictures below).

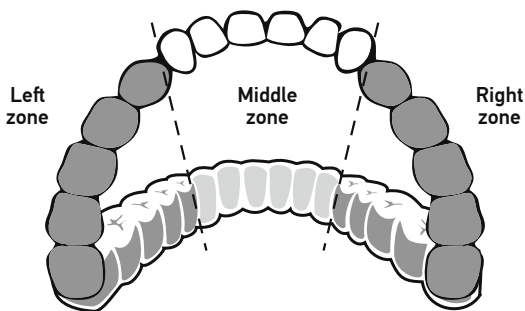
Repeat the brushing process for the lower jaw in reverse order (left/middle/right) as it is shown on the pictures.

#### Teeth brushing technique of left/right zones

#### Teeth brushing technique of middle zone

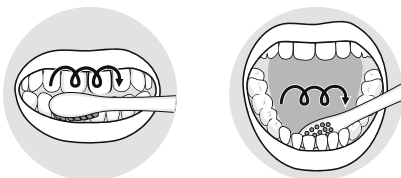Teeth brushing zones change every 20 seconds after a pause in the vibration of the brush.

## OPERATION VIA MOBILE APPLICATION

### ! NOTE

Mobile application provides opportunity to:

- teach your child to brush his/her teeth properly in a playful way together with the hero of the Brushy app;
- develop a healthy habit of regular tooth brushing;
- keep statistics on the quality of brushing and get recommendations on which teeth should be brushed more carefully.

**Install the application:** Android 5.0 and higher, iOS 10.0 and higher.

**Option a:** open the PlayMarket or AppStore, enter «AGU Brushy» in the search box and download the application.


**Option b:** for direct link to download the application, scan this QR code.

- 1 Download and run the application.
- 2 Register or log in through social media site.
- 3 Select device to connect.
- 4 Enter the required user data.
- 5 Switch on the device and press «Connect» button.
- 6 Follow the application instructions.

We start brushing teeth with the upper jaw on the right side.

During teeth brushing we go through three zones consistently (right/middle/left), each in accordance with it's own brushing technique (see pictures below).

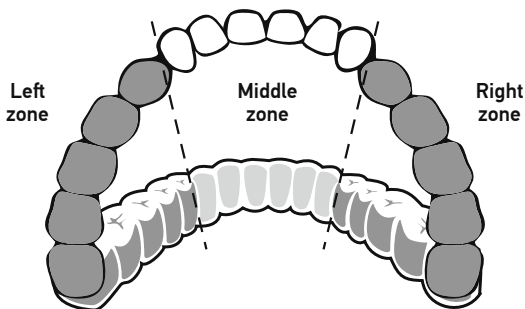
Repeat the brushing process for the lower jaw in reverse order (left/middle/right) as it is shown on the pictures.

## Teeth brushing technique of left/right zones

## Teeth brushing technique of middle zone

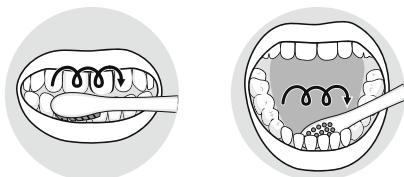
Teeth brushing zones change every 20 seconds after a pause in the vibration of the brush.

## 10 DEVICE CLEANING

Rinse the toothbrush head under a stream of water for a few seconds after each use.

If necessary, wipe all surfaces with a soft cloth. Do not use sharp objects for cleaning.
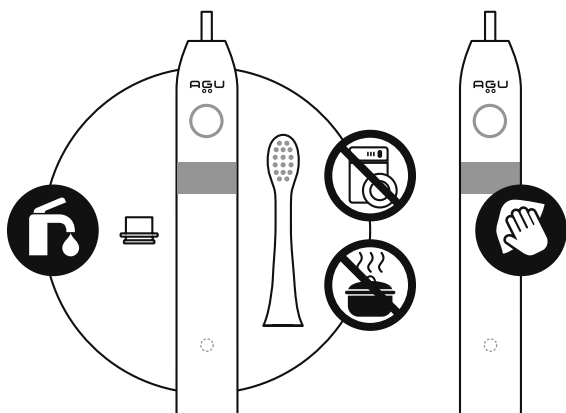
### **!** WARNING

Do not put the heads and the handle in the dishwasher.  
Do not wash the device in boiling water.

## 11 TROUBLESHOOTING

Description	Potential causes	Solutions
The device does not connect to the application	Internal battery low	Charge the device
	Wireless data transfer on the smartphone is turned off	Enable wireless data transfer on your smartphone
The device does not turn on	Internal battery low	Charge the device

## 12 STORAGE, TRANSPORTATION AND OPERATION RULES

- The device should be transported within the temperature range from  $-10^{\circ}\text{C}$  to  $+40^{\circ}\text{C}$  with relative humidity up to 75%.
- The device should be operated within the temperature range from  $+10^{\circ}\text{C}$  to  $+40^{\circ}\text{C}$  with relative humidity up to 65%.
- The device should be stored within the temperature range from  $+10^{\circ}\text{C}$  to  $+30^{\circ}\text{C}$  with relative humidity up to 65%.
- Do not expose the device to thermal shock.

### ATTENTION

After transportation or storage at low temperatures, it is necessary to keep the device at room temperature for at least 2 hours before switching on.

## 13 DISPOSAL

This symbol on the product or its package means that this product does not fall under the category of domestic waste.

- Proper disposal of the device will prevent adverse environmental and human health effects.
- In order to protect the environment, the device must not be disposed of together with domestic (household) waste. Disposal shall be provided in accordance with local regulations.
- The device must be disposed of in accordance with the EU Directive 2012/19/EU WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local public utility responsible for waste disposal.

## 14 CERTIFICATION

Smart kids toothbrush **AGU SKT6** complies with:

- the basic requirements of Radio Equipment Directive 2014/53/EU; Directive 2011/65/EU RoHS on the restriction of the use of certain hazardous substances in electrical and electronic equipment;
- the basic requirements of TR CU 020/2011 «On electromagnetic compatibility of technical devices»;
- the basic requirements of TR CU 007/2011 «On the safety of products for children and adolescents».

## 15 MANUFACTURER'S WARRANTY

The warranty period of this product is 24 months from the date of sale, providing that all operation conditions specified in this instruction manual are strictly observed. The warranty is valid only on presentation of the sales receipt and the warranty card completed by the authorized dealer. The warranty does not apply to consumables, constituent parts and components, such as battery, package, etc.

Warranty and free maintenance service is not provided in case of:

- violation of operating instructions;
- damage caused by intentional or erroneous actions of the consumer due to improper or negligent handling;
- evidence of mechanical impact, dents, cracks, chips, etc. on the device housing, evidence of opening the device, disassembly, unauthorized repair, ingress of moisture, exposure to aggressive agents or any other unauthorized alterations of the device, and in other cases of violation of storage, cleaning, transportation and operation rules specified in this instruction manual.
- ingress of oils, dust, insects, liquids and other foreign objects inside the device.

Follow the instructions precisely to ensure reliable and long-term operation of the device.

In case of abnormal operation of the device, please contact your retailer or Customer Service.

The manufacturer reserves the right to make alterations in the design of the device.

### **ATTENTION**

Do not try to repair the device yourself – this will void the warranty.

The device alterations without the manufacturer's permission will void the warranty.

For repair and maintenance, please contact a specialized after-sales service (see [agu-baby.com](http://agu-baby.com)).

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Montex Swiss AG,