


**AGU**<sup>®</sup> ADVANCED GROWING UP

**AGU**<sup>®</sup> ADVANCED GROWING UP

**AGU**<sup>®</sup> ADVANCED GROWING UP

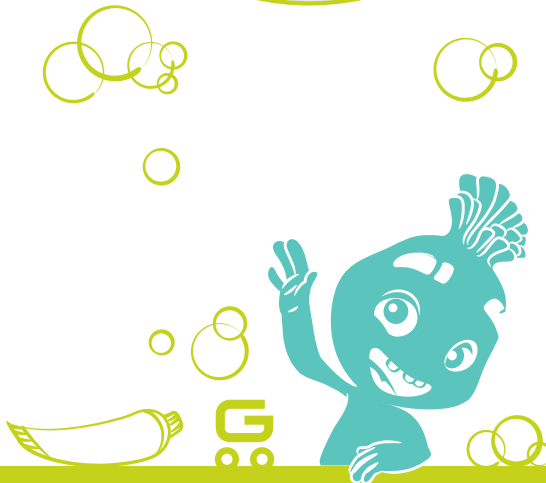

Brushy  
**SMART KIDS TOOTHBRUSH**


agu-baby.com

AGU<sup>®</sup> is the registered trademark by Montex Swiss AG,  
Tramstrasse 16, CH-9442, Berneck, Switzerland

agu-baby.com


swissdesign

# 1 PARTS DESCRIPTION

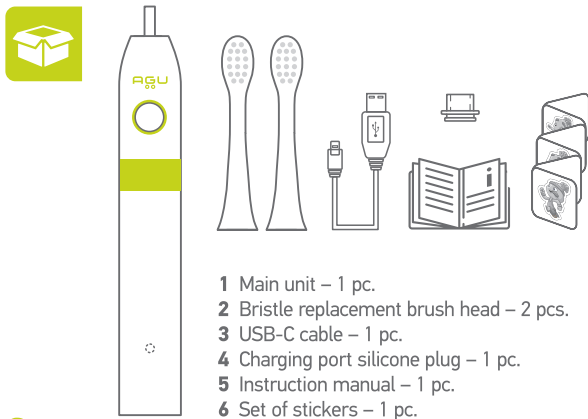


- 1 Main unit – 1 pc.
- 2 Bristle replacement brush head – 2 pcs.
- 3 USB-C cable – 1 pc.
- 4 Charging port silicone plug – 1 pc.
- 5 Instruction manual – 1 pc.
- 6 Set of stickers – 1 pc.


# 2 POWERING

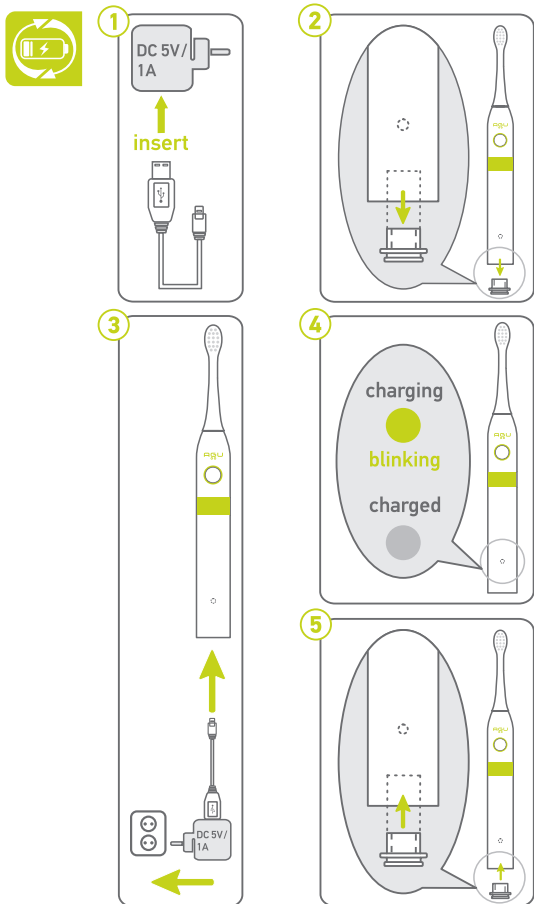

# 3 BRUSHING TECHNIQUE

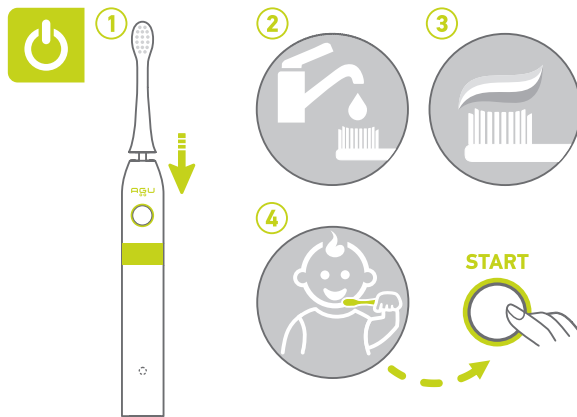

# 4 OPERATION OF THE DEVICE WITHOUT THE APP

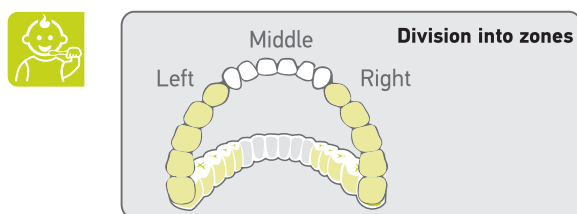


- 1 We start brushing teeth with the upper jaw on the right side.
- 2 During teeth brushing we go through three zones consistently (right/middle/left), each in accordance with it's own brushing technique (see pictures below).
- 3 Repeat the brushing process for the lower jaw in reverse order (left/middle/right).

Teeth brushing zones are changing every 20 seconds after pause in brush vibration.

## Teeth brushing technique, left/right zones



## Teeth brushing technique, middle zone

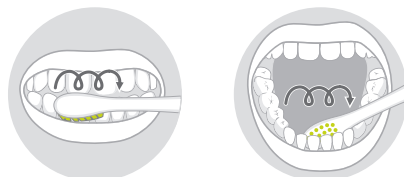

# 5 OPERATION OF THE DEVICE USING THE APP



# 6 MAINTENANCE



**Warning!**  
Do not wash replacement brush heads, main unit in a dishwasher.  
**Warning!**  
Do not wash the device in boiling water.

